



Surgical Back-Support-Corset according to Professor Lindemann  
(elastic 2/3 length corset, unisex)



# Surgical Back-Support-Corset according to Professor Lindemann

(elastic 2/3 length corset, unisex)



Compression of the abdominal cavity will show an improvement in posture, the posture will be more erect. The spinal column as well as the back muscles are relieved. Due to a reduction in pain the protective posture will be replaced by more physical activity as well.

## INDICATION (according to Prof. Petersen)

- static muscular weakness/osteoporosis, in conjunction with shoulder-straps in the area of the dorsal vertebra
- osteochondrosis, spondylosis, spondylarthrosis, mild scoliosis, basopobia, 1st. degree spondylolisthesis
- limited use after lumbar nucleotomy
- not recommended during inflammatory processes or 2nd. and 3rd degree spondylolisthesis.

## ORDER

Surgical back-support-corset according to Professor Lindemann or likewise elastic 2/3 length corset.

## COLOURS / SPECIAL FEATURES

The surgical back-support-corsets can be ordered in the following fabrics:

- 1.) Twill, 65% Dacron, 35% cotton, skin tone
- 2.) Jacquard, 100% cotton, skin tone

The lining is 100% cotton. All fabrics can be laundered at below 30° C (85° F)

## CLOSURES

Standard closures are: Hook and eye on the left, lacing front or side.

A more convenient alternative is the Velcro® closure, especially for patients of advanced age who may have difficulties with their hands.

We will of course manufacturer closures to your specifications.

## SPECIAL OPTIONAL EXTRAS



Inner-belt made of 20cm (8") wide lined tulle elastic mesh, has hook and eye closure.



Bottom closure (jock), can be installed later and is removable by way of Velcro® fastener.



Shoulder-brace with padded shoulder straps.



Support-strap with Velcro® fastener.




straps with Velcro® fastener

To finish the corset we offer the following:

Elastic edging, trimming ribbon in 21 and 25 mm width (4/5 and 1"), garters (narrow/wide), flexible spiral bones and surgical steels.

#### FINISHING TIPS

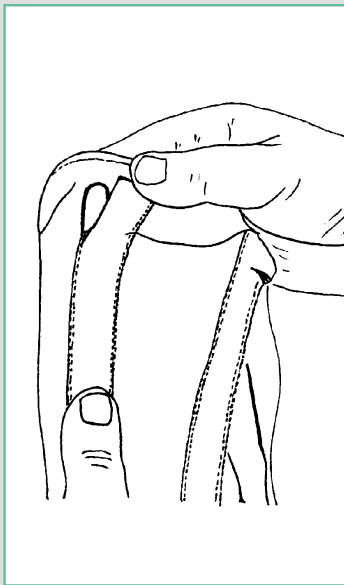

To facilitate laundering of the corset, we recommend making the surgical steels removable by leaving the top portion of the steel's pocket open and covering the open pockets with a flap.

Note:

Your surgical back-support-corsets will receive a personal touch by attaching your company's label.

Last but not least: We will realize your ideas and of course work to your specifications.

Give us a try!

**Figesta**  
CORSETRY MANUFACTURER

Ludwig-Erhard-Strasse 19 • 30982 Pattensen  
Tel.: +49 (0) 51 01 / 99 96 - 0 • Fax: +49 (0) 51 01 / 99 96 - 66